



# FUTURE CARE LOW FREQUENCY

“THE END OF ALL  
DISEASES”

## Contents

- Precautions before use
- Design and Concept
- Function
- Abdominal Doza
- Foot Massager
- Product composition
- Product parameters.
- Operation interface instruction and usage
- Remote controller and socket instruction
- Product connection
- Usage method
- Product usage
- Notes for use
- Fault diagnosis and treatment method
- Maintenance and cleaning

## Precautions before use

- Make sure to follow the product instruction manual before opening or using.
- Confirm whether the voltage is consistent.
- Use the power plug that meet the specifications.
- Do not connect the main power supply when the external equipment is installed on the body.

## Design and Concept

- It is designed as a high grade product and is easily detached due to its smart structure.
- It is ergonomically designed for the convenience of users.
- It is designed to image a soft human body of “s”line shape when viewed from the
- behind or the sides of the product.
- It is designed to provide diversity of choice.

## Function

- It is detachable and compact.
- It is high grade article and light.
- It provides maximum benefits as an excellent programs can be chosen in each mode.
- It is possible to upgrade programs and to install many other programs.
- Special feeling and special waveform.
- It is possible for 4 people to use the product at the same time.
- It has enhanced functions for the foot and each body part.
- large numbers of programs are available in each system with patented Doza.

ed from the



## Abdominal Doza

- It is ergonomically designed to increase usage effects.
- It increases treatment efficiency by computerizing the effect of therapy by each body part.
- It contains new materials to increase therapeutic efficacy and enhance deodorant function.
- It has patented system circuits inside to install many programs.



## Foot Massager

- It is detachable for increased usability.
- It has an external structure for less A/S.
- It increases wave transmission force and acupuncture with the standing type therapy.
- It maximizes treatment effects by using heating.
- It increases the foot correction effect by excellent adhesion regardless of different foot sizes.
- It increases usage effects with various output patterns.



## Product composition



Support master

Operation master

### Basic parts



Remote control

Abdominal Doza

Footing

handles

acupuncture patch


## Product parameters

Name of product	Health Frequency Therapy Device
Name of model	FUTURECARE-7
Power supply	AC100-240V/50-60Hz/60W
Output wave	Biowave
Maximum current	10mA(Auto-Control)
Size of the main body	330(L)×250(W)×60m(H)
Size of the footing	350(L)×310(W)×35m(H)
Main Weight	2.1 Kg
Whole set weight	9Kg




# Operation interface instruction and usage

# Remote controller and socket instruction

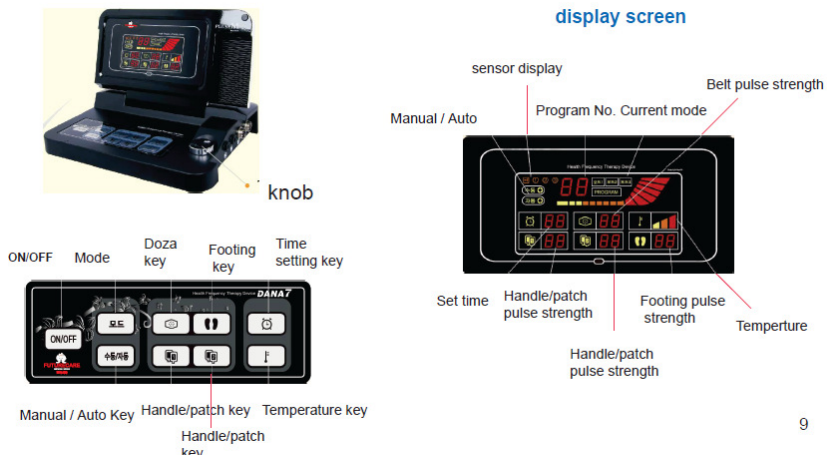

9

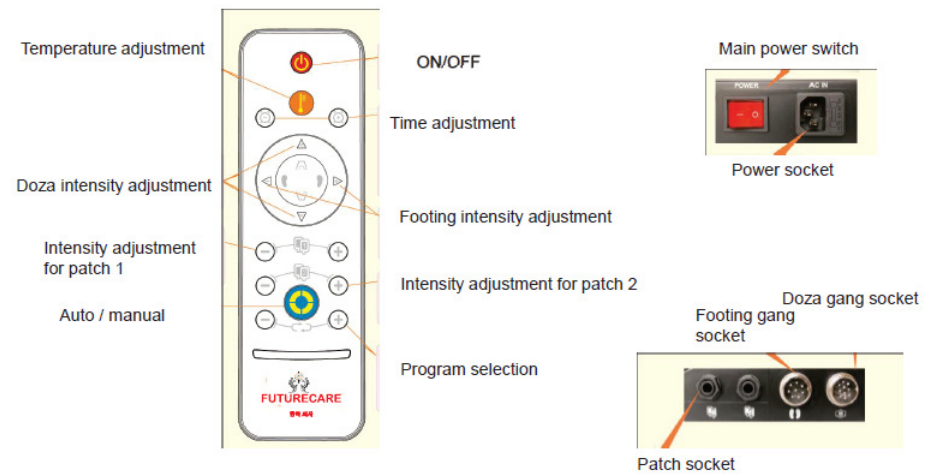

## 6 // Product connection

### ○ Basic connection

#### 1. Power connection

Connect the power cord to the management instrument, then connect the other pin to the socket



#### 2. Doza connection

Inserted the Doza correctly in the socket groove, as shown in the attached picture and then tighten the screws



#### 3. Footing connection

The foot connecting groove is upward, and the connecting wire is connected with the socket on the foot connection as shown in the figure



#### 4. Pad connection

After the two lines of the pad are connected on the respective pads, the back bolt is connected with the socket on the pad



### ○ Optional accessories connection

#### 1. Cushion connection

The connection socket can be connected with the cushion, and the groove is upward



#### 2. Pads and handles connection

The pad socket can connect pad or handles as shown in the picture, you can choose one to use



## 7 // Usage method

### ○ Doza

Check the Doza connection socket is connected with the doza, apply the electric gel to the carbon part inside of the doza. Then tighten the doza and fix its management position.



### ○ Footing

Check the footing connection socket is connected with the footing, apply the electric gel to the carbon part inside of the footing. Adjust the height of the chair and the distance on foot completely close to the heel part of the footing, relax the feet and body.



### ○ Pad

Make sure that pad is connected to the correct socket, tighten the belt, fix the treatment area, not too tight or too loose




### ○ Cushion and Pad

Confirm the connection state of selected components and connecting sockets, apply the electric gel to the corresponding carbon part. Put the cushion on the floor, make the body close to it. The pad is fixed tightly by the belt.



○ Patch usage for different body part

Generally, the electrodes are arranged at the treatment site. Please consult your doctor if putting the patch near heart

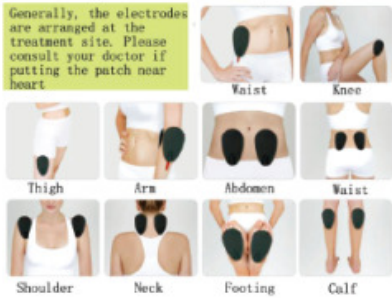

○ Doza usage for different body part



8. Product usage

○ Intensity adjustment for doza and patch



## Product usage

- The main power from OFF state changes into ON state.
- Touch the operating ON/OFF button ⑦ to become ON state.
- Make sure that all output window(A-D) levels are 00 state.
- Touch the button with a belt pattern on it, adjust the intensity to 1-20 by knob button. Then the window (A) display output intensity.
- Touch the button with a foots pattern on it, adjust the intensity to 1-20 by knob button. Then the window (B) display output intensity.
- Touch the button with a pad1 pattern on it, adjust the intensity to 1-20 by knob button. Then the window (C) display output intensity.
- Touch the button with a pad2 pattern on it, adjust the intensity to 1-20 by knob button. Then the window (D) display output intensity.
- Touch the button with a clock pattern on it, adjust the intensity to 1-20 by knob button. Then the window (D) display remaining time.
- Touch the button with a thermometer pattern on it, chose H, M, L grade. Then the window (F) display corresponding temperature strength.

## Product usage

### Time and temperature adjustment

The temperature should be adjusted to an appropriate temperature before use.

The temperature can be adjusted to H, M and L three grade.

Time can be adjusted to 1-90 minutes, and can be extended or shorten according to the using time.

It will give an alarm once the set time come, and automatically cut off the output. All of the output intensity change to 0.

### Knob usage

If you press each function button, output light flashes, then you can rotate the knob to set the output intensity.

Stop rotating the knob after about 5 seconds, the output light stop flashing.

### Mode adjustment

Different mode can be set by selecting program

Program1:1-10 mode.

Program2:11-20 mode.

Program3:21-25 mode.

### Auto / manual setting

Auto / manual buttons can be used to set automatic or manual.

Manual: can use the setting program continuously.

Auto:Automatic convert program according to the mode.

For example:When selecting program 1, 1-10 mode cycles automatically.

## After using

- After using,must press the ON/OFF key and turn off the main power switch. Wipe the gel of the Doza and footing by wet towel.
- This is the best way to maintain the state of the product, please operate as required.

## After-Sales Service

After-sales service must go to the shop where to buy the product or service center to consult.

Disassembly can not be arbitrarily disassembled products, if the non professional disassembly of any decomposition of the product, the company will not be responsible for the product after-sales service,Please go to the service center to consult or repair.

## Caution:

This product can be automatically converted to 90V-260V. However, when the voltage is changed, the fuse installed on the power socket must be replaced with the standard voltage.


## Fuse replacement method

- Remove power plug.
- The power supply socket and fuse are integrated. The fuse is installed in the lower end of the power supply socket.
- Open the fuse box by a straight screwdriver.
- Replace the fuse as the picture below, let the power socket highlight the part up and push it in.



## Remote controller battery replacement method

- ① Push the cover on the lower end of the Remote controller according to the direction of the arrow, and detach it.
  - ② Install the battery according to the positive and negative electrode.
- ※ This product uses 2 AA batteries.

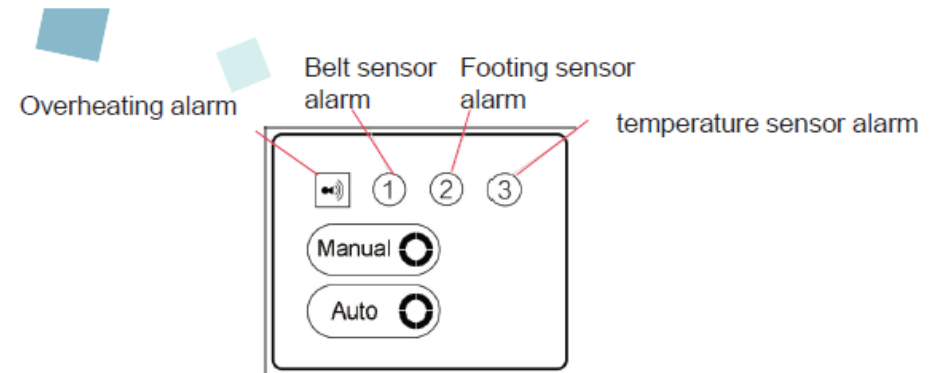

- Belt sensor alarm: the alarm light is off when connector is normally connected to the belt. If alarm light is on when the belt connector is properly connected, please consult the after-sales service center.
- Footing sensor alarm: the alarm light is off when connector is normally connected to the footing. If alarm light is on when the foot connector is properly connected, please consult the after-sales service center.
- Overheating alarm and temperature sensor alarm light: If overheating alarm and temperature sensor alarm light on, it shows the machine already cut off the internal current to protect the user from the risk of burns by the belt or the foot of the high temperature.

## Safety precautions

Must read "safety precautions" before using  
The instructions record the important content of safety, if not followed the instructions, there may be several dangerous degree happen.



**If you ignore this sign and do the wrong operation, there will be a possibility of injury or damage the machine.**

✘ Using the product correctly and safely is to take preventive measures. Please follow the rules.



## Prohibition of use

1. Persons who have the following circumstances shall stop using:  
Patients with acute disease.  
People with high fever.  
Pace Maker.  
Perceptual obstacle.  
People in pregnancy.  
Patients with malignant tumor.  
Operation  
People who use artificial organs and artificial bones.
2. Prohibit the use of the following body parts  
Eyes.  
Genital organs.  
heart.  
brain.  
vocal cord malignant tumor.  
Plastic surgery part.  
Wound around area, head and ears.  
Avoid contact with teeth directly.  
The part where insert metal, plastic and silicone parts.
3. Patients in the hospital to consult their doctors.



## About product using

- There will be a risk of Machine damage or abnormal if decompose, repair or transform the machine arbitrarily
- Stop using the product immediately if product abnormal, consult the service center or the place of purchase.
- Use from low to high intensity (do not use high intensity suddenly, the body will not adjust).
- Adjust the intensity to the level of the user's skin.
- Not available for purposes other than the use of this specification.
- Stop using the product immediately if you feel uncomfortable, consult your doctor.
- Do not throw the product, or give a strong impact on the product.
- Do not grab the cable and rotate the belt.
- Do not use the product on the eyes, it's very dangerous.
- Do not use the product with other devices.
- Do not use in the bath or damp place.
- Use the product no more than 3 times.
- Use 3 times a day, the interval must be maintained at more than 4 hours.
- Do not go out or move the machine to another place when the machine working.



## About power

- This product is AC90V-260V.
- The power plug must use a socket with Earthing-contact, there will be a risk of electric shock when fault or leakage.
- Do not plug a plug on multiple sockets. There will be a risk of fire when socket heating.
- Do not pull the power plug (There will be a risk of fire when bolt damage).
- Do not damage, bend, pull, tight or bind the power socket. It will cause fire or electric shock when socket or cable damaged. Please use the cables, wires and bolts properly.
- The power socket must be kept clean and no foreign matter entered (There will be a risk of fire when entered foreign matter).
- There will be a risk of electric shock or abnormal if decompose, repair or transform the machine arbitrarily. (If the product is transformed privately, even if within warranty period, the company is not responsible for free Maintenance).
- Do not use wet hand to plug the power supply or unplug the power, there will be a risk of electric shock or fire.
- Do not use damaged socket or bolt, there will be a risk of electric shock or fire.
- If the cord is damaged, it must be replaced by the manufacturer and its maintenance department in order to avoid danger.
- Do not put the power socket on the walking channel.
- Unplug the power cable when it thunder.
- Unplug the power cable when not use.



### About product installation

- Don't knock the product or give the impact when handling the product (May damage the product).
- Avoid direct sunlight (May change the original color).
- Do not close to flammable products such as alcohol (May lead to deformation or fire).
- Do not install the product on the inclined section or shock location (installed it on a flat place).



### About children:

- Play around the product may encounter the equipment and get wounded.
- The children are not allowed to close when using the product.(if not noticed, the child will be frightened when the machine starts).
- Do not allow children to play with this product (there will be a risk of suffocation if cord around the neck).
- There will be a risk of electric shock or fire if bite the wire or scrap with a knife.
- Cut off the electricity supply when not use (children may be frightened when touch it).



### About product usage

- Must be familiar with the usage and methods of operation.
- Don't touch sharp objects.
- Placing beverage around the product may tilt and affect the device starting, even lead to electric shock accident.
- If the belt or the instrument be pured into the beverage or water, etc. unplug the power immediately and contact after-sales service center.
- Do not put pressure on the product, it may be damaged .
- Product broken or power failure, unplug the power immediately and stop using.

## Fault diagnosis and treatment method

Before repairing, please check with below methods. Please consult the after-sales service center once confirm abnormal operation.

Symptom	Inspecting Point	Treatment method
Machine not start	<ul style="list-style-type: none"> <li>•Whether the power switch on the socket is off</li> <li>•Whether the main power on</li> <li>•Whether the ON/OFF switch on</li> </ul>	<ul style="list-style-type: none"> <li>•Confirm whether the voltage is consistent</li> <li>•Plug in the socket</li> <li>•Turn on the main power</li> </ul>
Remote contraller not work	<ul style="list-style-type: none"> <li>•Whether the batteries used for a long time</li> </ul>	<ul style="list-style-type: none"> <li>•Replace to new batteries.</li> </ul>
The intensity is too low	<ul style="list-style-type: none"> <li>•Different parts, the output intensity is also different</li> <li>•Whether to turn outside in the pad to use</li> <li>•Whether to connect the connector of the using part.</li> <li>•Whether adjust the output intensity properly</li> </ul>	<ul style="list-style-type: none"> <li>•Relate to the physical condition of individual</li> <li>•Relate to different seasons</li> <li>•Make the body close to the treatment area</li> <li>•Connect the connector of the using part correctly</li> <li>•Adjust the output intensity properly</li> </ul>
Feel stimulated when using	<ul style="list-style-type: none"> <li>•Whether adjust the output intensity properly</li> </ul>	<ul style="list-style-type: none"> <li>•This product is a soft feeling product, if you feel uncomfortable, please stop using immediately.</li> <li>•Apply the electric gel to the using part before using.</li> <li>•Please consult the after-sales service center if still not work after doing the treatment above</li> </ul>
Doza and footing not heat	<ul style="list-style-type: none"> <li>•Whether open the temperature adjustment button</li> </ul>	<ul style="list-style-type: none"> <li>•Adjust the temperature adjustment button properly</li> </ul>



## Maintenance and cleaning

- After using, store the wires, belts, and other components separately to avoid loss.
- Avoid direct sunlight and prevent dust into the product.
- Stored at cool, dry and ventilated place.
- Don't knock the product or give the impact on it.
- If the product is transformed privately, even if within warranty period, the company is not responsible for free maintenance.
- Products should be kept clean after use.
- Prevent tilt and vibration. maintain the product in a safe state.
- After a long time not to use, please comply with the rules when start again.
- Unplug the power before cleaning the machine.
- Use a soft cloth to wipe the body (especially the display screen). don't wash with water, do not use rough cloth.
- Wipe the dust off the mat with a soft and dry cloth.
- Do not use benzene, thinner or alcohol to wipe the machine and basic parts.
- Fully dry after cleaning.

## Warranty \_\_\_\_\_

We give you the following limited warranty for this product.

<b>Product Name</b>	<b>Low Frequency Stimulator</b>
<b>Model Name</b>	<b>FC-7</b>
<b>Customer</b>	<b>Address:</b>
	<b>Name:</b>
	<b>Mobile:</b>
<b>Dealer</b>	
<b>Date of Purchase</b>	
<b>Duration of Warranty</b>	<b>Within one year from the Date of Purchase</b>

- This is the product which has passed strict quality control.
- This warranty should be well retained.
- If you have any question during operation, please get in touch with the dealer you bought the product from.
- We will provide you with a kind reply in detail.

This limited warranty promises you to apply to the malfunction occurred during the period from the date of purchase and only the part warranted as described in Warranty regulations as under. Please refer to <Warranty Regulations> for further information.

## Warranty Regulations

1. In case the product is out of order while working under normal conditions as per the Operation Manual, the Warranty will apply to repairing within 2 years from the date of purchase with free of charge. It, however, will be charged after the warranty period.
2. As cords (cables) and other parts are the consumptions, warranty period is within 6 months from the date of purchase. After the period, it should be paid.
3. If the product is out of order within the warranty period, the product and this given warranty will have to be presented.
4. The given Warranty should be filled up, unless otherwise it will be invalid. Make sure to fill it up at the time you buy the product at dealer shop.
5. The given Warranty will never be reissued, the product will not be repaired if the Warranty is lost. The Warranty, therefore, should be well retained.
6. Following are the objects to be charged even though it is within warranty period:
  - a. Malfunction by improper repairing or remodeling.
  - b. Malfunction by Shocks from exterior like dropping.
  - c. Malfunction by infiltration of liquid or dust.
  - d. Malfunction occurred when in use under circumstance against Warranty Regulations.
  - e. Malfunction caused by natural calamity like fire, earthquake, flood and thunderbolt, or pollution or abnormal electric voltage.
  - f. In case the given Warranty is not Presented.
  - g. In case date of purchase, customer name and dealer name are not filled up or are changed in the given Warranty.