



FUTURECARE

MEDICAL KOREA ○


한국 의사

CHILBOSUK CERAMIC HEALING BELT

Make sure to read this manual for the safety of customer and proper use of the product.

This manual includes the product warranty certificate.

The specification and design may be changed to improve the performance without a prior notice.

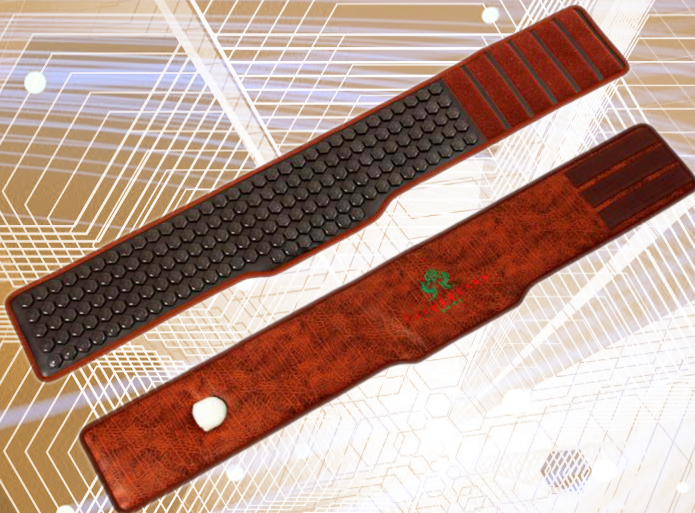

FDA Certificates

Welcome To FUTURECARE Family

**We Sincerely appreciate for our customers
who have purchased our products.**

We at FUTURECARE, have devoted our utmost effort in marketing high quality health related product for out customer. Since our research and development from 1998 onwards.

We have continuously upgraded our technical system in order to develop into a company at forefront for a health conscious culture. As a result, FUTURECARE proudly maintains senior position with a major market share for the thermal massage industry having presence in many states of our respected country INDIA. We are importing and manufacturing best quality medical device and thermal massagers for our customers. We are mainly importing products from Asian & European countries. All staff members and employee, including myself, keep in mind that this kind of achievement has been made possible only through our customers continuous love and support. We promise to dedicate our best efforts toward maintaining customer satisfaction as one of our top priorities. 21st century is rapidly emerging as an era that emphasizes the value of living a healthy life. Now, FUTURECARE will do its utmost in becoming a milestone for this changing lifestyle by accelerating product development and improving customer service based on our company's philosophy and pursuit better health. We are a company suitable for anyone and everyone. Please keep and eye on the changing aspects of FUTURECARE progress.

Thank you
Dr Jeon Hun Baek
CEO

FAR INFRARED CURES, VITALISES AND WORKS PREVENTIVE

What is far infrared and is it safe?

Infrared light or radiations are electromagnetic waves that form a part of the light spectrum. Sunlight contains 4% ultraviolet radiation, 52% infrared radiation and 44% visible light. The physical characteristics of FIR (Far infrared Radiation) are: Invisibility, Linearity, refractivity and reflectivity. Upon being absorbed by various objects, it has the ability to produce warmth. Far infrared harmonizes with the natural life-force released by human body to enhance resonance absorption and has the ability to penetrate deep into the human body. Far Infrared Rays are the healthy rays contained in sunlight, it's most essential and beneficial for life and growth of all living things on this planet. This is the reason why Far-infrared is referred to as, "The Light of Life". Far infrared is healthy and save, this stands in contrast to gamma, microwave or ultraviolet radiation. A large dose of this types of radiation could be severely harmful. It is the infrared heat in the sun that gives us a sense of comfort. It is proven in many studies that far-infrared rays support the growth and health of living cells in plants, animals and humans.

Far-infrared in our body

The human body also produces far infrared heat, as the result of food combustion. The intensity of our body heat fluctuates continuously. When food combustion reaches optimum level we feel healthy, being able to overcome ailments, viruses and bacteria. When food combustion does not function properly, it has a negative effect on cell metabolism and thus our general health. As a result, there is less released far infrared heat, which is important in order to keep the body warm and protect it's vital organs. We become more susceptible to disease and we tend to age faster. An external source of infrared heat can also be a major source to stimulate digestion processes. Of course, the body does needs enough (easily absorbed) food/power/fuel to generate healthy heat*.

** In case of illness and too little food intake, ultimately our (fat) reserves are addressed and if that's burnt-out then other body tissues will be used or exposed to oxidation. In other words, the body becomes exhausted and is getting sicker.*

FAR-INFRARED CURES, VITALISES AND WORKS PREVENTIVE

Far-infrared and Healing power

Our hand palms radiate far infrared energy from between 8 and 14 microns. In principle, it is possible to alleviate pain and cure ailments by means of hands-on healing. Why we move our hands automatically to that place where we are in pain? For some people the personal emission of FIR is developed by practicing Reiki, Pranic Healing, Yoga, Qigong etc. and they will be able to heal by laying on hands. No wonder... this form of medicine is already been used for thousands of years. Mother Nature has given us every opportunity to live a healthy life. Regular exposure to the far infrared heat from the sunlight makes us feel healthier and more energetic. We are nourished and our immune system sustains strong. It has been shown that the far-infrared waves between 6-14 microns have special regenerative effects on the human body. If a cell or tissue is exposed to these waves, then the reaction mechanism is optimized as well as the adhesion of nutrients and an optimized osmosis of water molecules through the cell membrane takes place. The immediate consequence is an increased cell membrane potential, an improved transport of photon particles and ultimately an improved cell division. These waves will take care of that the body's healing processes are stimulated optimally.

Far-infrared and the detoxification process

Due to today's high levels of stress and toxicity the body's natural ability to heal, is often compromised. Some of us have reached critical stress and/or toxic levels and we do not have the energy to detoxify or repair our bodies. Far Infrared Thermal Treatment can leave you feeling relaxed and energized. When toxic gases such as sulfur dioxide and carbon dioxide or toxic substances such as lead, mercury, or chlorine reaching the big water molecules, they are encapsulated. Where these toxins accumulate, blood circulation is blocked and energy supply of the cell becomes interfered. The detoxification process begins when 8-10 micron far infrared waves reach the large water molecules, heats them up and initialize them to resonate. These vibrations break the bonds of the ion-toxins with the water molecules. Now a division of large water molecule clusters ensues and releases the encapsulated gases as well as other toxic substances. Japanese scientists report that a similar process occurs with the use of far infrared therapy in case of clogged capillaries: the far-infrared heat dilates the fine capillaries where the process starts to dissolve trapped toxins and waste.



FAR INFRARED CURES, VITALISES AND WORKS PREVENTIVE

Chilbosuk®

The **Chilbosuk®** stones are baked ceramic stones containing Tourmaline crystals, Geranium, volcanic stone and other trace mineral. In 1880 it was discovered that Tourmaline has a very special crystal structure. It is the only stone that under the influence of solar power constantly emits a positive and negative charge at each end of its crystals. This phenomenon is called pyro-electricity (0,06A). It is a form of electricity that is healthy and increases metabolic activity at a

cellular level. Furthermore Tourmaline emanates solar infrared (FIR) radiation and negative ions, which intensify further during the heating process of the stone. Lying on the heated **CHILBOSUK®** stone has a beneficial effect on the respiratory system, improves circulation processes and revitalized the cells. It helps to keep the body healthy and promoted the discharge of carbon dioxide, heavy metals, and other waste products from the cells. The far-infrared rays of these stones are not only beneficial for the skin and muscle tissues, but for all cell tissues of the body and have also an anti-bacterial effect. Even with chronic problems and illnesses one can experience relieve of certain symptoms and gain an improved condition.

About Positive and Negative Ions

An ion is an electrically charged atom, molecule or other group of bonded atoms. An ion can be positively or negatively charged respectively by a lack or an excess of one or more electrons. The atmosphere around the earth is filled with positive and negative ions. The earth itself has a negative charge, so negative ions are repelled. This repulsion ensures that there is a certain balance between negative and positive ions in a ratio of 10 to 12, respectively. There are always more positive than negative ions in the average environment. Negative ions vitalize the air, making the air healthy and livable, thereby having a beneficial effect on all living things. Negative ions are therefore considered as vitamins of the air. An abundance of positive ions decates the quality of the air and is seen as a major cause of both physical and mental health issues, ranging from headaches to depression and even heart attacks. Research has shown that negative ion levels above 1000 ions per cc have beneficial effects on humans and other organisms. Unfortunately, negative ions end up into ever smaller concentrations because the conditions under which negative ions are produced are greatly reduced by human activities such as urban development and environmental pollution.



CHILBOSUK CERAMIC HEALING BELTS



Product Specifications

BRAND NAME	Futurecare
Model Name	Chilbosuk Ceramic Healing Belt
Model No.	
Manufacturer name and address	<Futurecare> 124 MYEONGJANG-RO 63 BEON-GILL, DONGNAE-GU, BUSAN. SOUTH KOREA
Product License Number	330100110194
Weight or packing unit	
Other necessary matters	1) Constant frequency and rated voltage of the rated power supply : 220v 50/60 Hz. 2) Power consumption at the maximum output : 90 W. 3) Type and degree of protection to the electric shock : Class 2 Equipment Type B Applied part. 4) Range of heating temperature : 30°C-70°C. 5) Operating time : automatic off after 8 hours of consecutive output.

Chilbosuk[®]

Temperature settings

35-40°C / 95-104°F | Safe for Sensitive + Extended Use | *Regenerative Sleep + recovery; Gives deep relaxation; Combats Jet Lag; Helps digestion; Headache and migraine reduction; Improve immunity; Improves energy levels; Regeneration of the cells.*

45-50°C / 113-122°F | Cycles of 2 Hours or More | *Sports + Flexibility Training; Cardio + Respiratory Aid; Blood Pressure + Sugar Regulation; Improves glandular functions; Female Cycle Balancing + Support; Improves hormonal balance; Use with asthma, allergy symptoms and respiratory problems.*

55-60°C / 131-140°F | Cycles of 30-90 minutes | *Helps with Stress, Tension, Anxiety + Trauma Relief; Increases resistance to stress; Joint, muscle and weak tissue support + Pain Relief; Profound Relaxation + Mood Elevator; Increases vitality; Improve immunity + Detox.*

65-70°C / 149-159°F | Cycles of 30-60 minutes | *Sauna Level Cleansing + Detoxification of heavy metals and cellular waste; Improve viral, cystic, acidic and endemic malfunctions; Natural Weight Loss + Skin Tone; Enhanced Sexual Vitality + Metabolic Aid in case of low immunity, flu, cold, hangover and low energy levels; Healthy Blood Circulation + Purification; Environment + Bedding Sanitizer.*

And further:

- Heat up the mattress or pad 10 minutes before use.
- Far Infrared can help the self-healing ability of the body to become reactivated again.
- It is of great importance to drink a minimum of 2 liters of additional water in a day (no cold drinking water, drink at room or body temperature).
- Use thin cotton clothes or put a thin cotton sheet over the mat. Rather no synthetic clothing (while sweating).
- In case of pregnancy and for people with a pacemaker apply not above 42°C / 108°F and no more than 1 to 1½ hours use afterwards.
- Do not use in fast growing malignant tumors or serious skin injury.
- Check with specific complaints always your physician/practitioner or therapist.
- Children and elderly people should use our products always under supervision.
- Keep in mind that; (chronic) disease, ailments and limitations usually arise in a day and so also not just improve or cure. Our products are not medicines health articles.

Greetings

Thank you for purchasing Chilbosuk Ceramic Healing Belt. We believe that this product will help you lead a healthier life. To use the products safely for a long time, please read this manual thoroughly and follow instructions on how to operate the product.

Exemptions

The manufacturer, importer and seller are not responsible for any problem or damages caused by user's in cautions or maintenance of the product. The manual is based on the current specifications of the model, but the product purchased may have different specifications since Futurecare Korea, is continuously adding new features and applying new technologies dimensions and functions may be changed to improve the performance without a prior notice.

Caution for Safety

Description

Caution !

In the event that there is a possibility of inflicting minor injury or property damage if not complied with the contents

Risk



In the event that there is a high possibility of inflicting serious injury or death if not complied with the contents.

Warning !

In the event that there is a possibility of inflicting injury or death if not complied with the contents

Prohibition



In the event that there is a serious possibility of inflicting death if not complied with the contents absolutely prohibited

Risk and Prohibition

1. Avoid pulling out the power cord or touching the power plug with the wet hand.
2. When there is thunder, lightning, power breakdown, or earthquake, stop using the product and disconnect the power immediately.(There is a risk of electric shock and fire.)
3. Absolutely do not use the product by women in late stage of pregnancy, spine ailment patient, patient wearing pace maker for heart beat, patients with severe heart ailment.
4. Do not use the product in closed space with fuel or chemical substance.
(It may cause the personal injury or fire.)
5. Do not use it in a place with substantial dust or severe change of temperature.
(It may cause the fire or electric shock.)

CONTROLLER FUNCTIONS



Temperature Display

Lamp comes on when you touch the dial.

Power Button

Temperature Control Dial

Turn the dial right and left to where the lamp goes off. Temperature may be set for maximum range, 30oC - 70oC (86oF - 158oF). Recommended Temperature is 50oC(122oF)




Coarse Rocks



Power



Fine Power

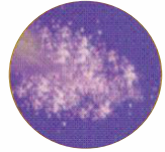


Careful and rigorous selection of high quality Korea produced Germanium and Elvan rocks, imported Volcanic rocks and tourmaline is made


Each of the four rock types is ground into fine power.

Tourmanium power is vacuum dried and further ground into fine power. It is then sifted by wind to collect only the finest micro-granular power.

Shaping

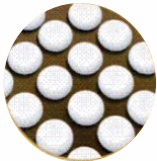

Baking

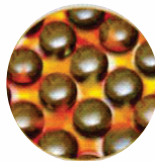

Chilbosuk Ceramic



Micro-granular tourmanium power is formed into molds for shaping.

Molded tourmanium is baked in over at over 1,300 C to remove any impurities.

The produced tourmanium therapeutic ceramic, the finest in the world, gives off deep shiny color without any coating and generates Far-infrared Ray Radiant Heat for therapeutic effects on



Benefits from Finger Pressure & Moxibustion
 Futurecare ceramic mat (personal electric potential stimulator by warm heat) has a multitude of round shape Tourmanium ceramic disks arranged densely in smooth contour so that they provide the optimum finger pressure and Moxibustion stimulation to your body during the product use, thus product used. Thus same effect as performing Moxibustion on the 365 meridians in

Far-infrared Ray Warm Heat Therapy
 Since the early years of Korea's long history, the Korean people have wisely employed the benefits of warm heat, which is Far-infrared rays that are emitted from clay walls, ceramic tiles, and 'Ondol' heated floors. Futurecare stimulator used this very principle so that the Far-infrared rays produced by its Tourmanium Ceramic array penetrated deep into your body to enhance metabolism and blood circulation. The product is also effective in relieving muscle tension and pains.

Far-Infrared ray (for better blood circulation)
 Far-Infrared ray provides body vitalizing energy. The radiant heat in the form of far-infrared ray generated in the Tourmanium ceramic penetrates deep into our body to enhance body metabolism and to stimulate blood circulation so that blood reaches to the finest blood vessels, giving your body vitality.

Anionization Function
 when the "Qi" (life-giving energy) in our body loses its balance, our body, condition deteriorates rapidly. Tourmanium ceramic has anionization function built-in so that it provides vitality to modern humans whose bodies often lose Qi balance when exposed to many forms of pollution daily causing positive ion presence.

- Six functions of far-infrared rays on human body**
- Warm heat function: maintains body temperature at right level
 - Maturation function: nurtures body growth and maturation
 - Balancing function: helps maintain nutritional balance in the body
 - Moisture control: maintains body moisture level at a right level
 - Neutralization function: helps discharge of body wastes and neutralize body
 - Resonance function: stimulates decomposition of nutrients for nutritional balance and enhanced body metabolism

Anion Therapy
 anion therapy was first developed by Dr. takada of oriental university of Japan. Drs. Nehru and Jackman of Germany, Nobel laureates in medicine, proved the existence of anion paths that connect cells with other cells. The effects of Anions on the human body are cell revival, blood purification, coordination of the sympathetic nerve system, and immunity enhancement. Futurecare has this beneficial anionization function. The negative potential stimulator Rod (provided separately) has been designed to apply directly to part of body that experiences pain.

PRODUCT INTRODUCTIONS




FUTURECARE CHILBOSUK CERAMIC HEALING BELT

Components



HEATING CONTROLLER


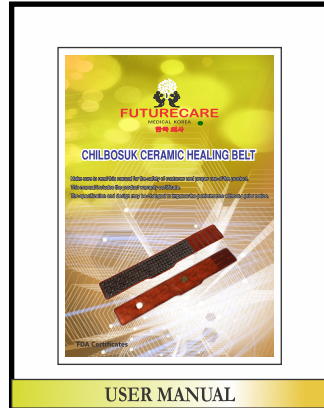


USER MANUAL



CHILBOSUK HEALING BELT

Uses of Healing Belt in Different Body Parts

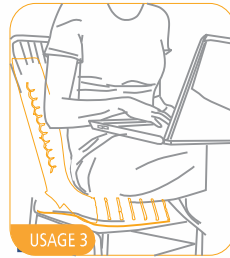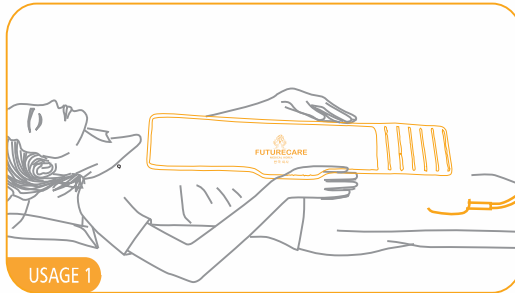

Product Introductions

Use of Chilbosuk Ceramic Healing Belt

Features

- Portable & Easy to Carry.
- Normalization of Intestinal Functions.
- Tourgia that emits Far-Infrared rays.
- Removal of Excess fat in lower Abdomen, Back and other parts of body increase metabolism and improve immune system.
- Relieves all body pains & tendos.

Note

- Use the Extension Belt of Futurecare Chilbosuk Ceramic Healing Belt is smaller than your size.

Warranty

We Gives You The Following Limited Warranty For This Product

Product Name	Futurecare Healing Belt
Model Name	FCB
Customer	Address:
	Name:
	Mobile:
Dealer	
Date Of Purchase	
Duration Of Warranty	Within One Year From The Date Of Purchase

This is the product which has passed strict quality control.

This warranty should be well retained.

If you have any question during operation, please get in touch the dealer you bought the product from.

We will provide you with kind reply in detail.

This limited warranty promises you to apply to the malfunction occurred during the period from the date of purchase and only the part warranted as described in warranty regulations as under. please refer to <warranty regulations> for further information.

Warranty Regulations

1. In case the product is out of order while working under normal conditions as per the Operation Manual, the warranty will apply to repairing within 2 year from the date of purchase with free of charge. It, however, will be charged after the warranty period.
2. As cords (cables) and other parts are the consumption, warranty period is within 6 month from the date of purchase. After the period, it should be paid.
3. If the product is out of order within the warranty period, the product and this given warranty will have to be presented.
4. The given warranty should be filled up, unless otherwise it will be invalid. Make sure to fill it up at the time you buy the product at dealer shop.
5. The given warranty will never be reissued, the product will not be repaired if the warranty is lost. The warranty, therefore, should be well retained.
6. Following are the objects to be charged even though it is within warranty period:
 - a. Malfunction by improper repairing or remodeling.
 - b. Malfunction by shock from exterior like dropping.
 - c. Malfunction by infiltration of liquid or dust.
 - d. Malfunction occurred when in use under circumstance against warranty regulations.
 - e. Malfunction caused by natural calamity like fire, earthquake, flood and thunderbolt, or pollution or abnormal electric voltage.
 - f. In case the given warranty is not presented.
 - g. In case date of purchase, customer name and dealer name are not filled up or are changed in the given warranty.

Futurecare Healing Belt



FUTURECARE

MEDICAL KOREA

한국 의사

KOREA OFFICE : 124 MYEONGJANG-
RO 63 BEON-GILL, DONGNAE-GU,
BUSAN. SOUTH KOREA, KOREA